

February 1, 2008

2D electrophoretic and other high-resolution separations and analysis of proteins in biological samples

**Helen Kim
934-3880
helenkim@uab.edu
McCallum Building, room 460**

Feb 1, 2008

HelenKim/UAB/PharmTox

1

Learning objectives

2-D gel proteomics: What it involves:

2D separation

image, statistical analysis

MS

Other types of 2D electrophoretic protein separations;

Free-flow



2D Blue native electrophoresis

Feb 1, 2008

HelenKim/UAB/PharmTox

2


Initially: The objective of protein separation in proteomics was to get a protein “spot” or “band”, for mass spectrometry analysis, to identify the protein/its modifications

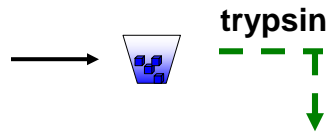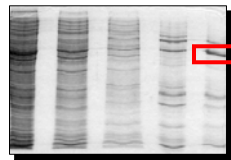

MALDI-TOF mass spectrometry

ID of parent polypeptide

MASCOT search engine
<http://www.matrixscience.com>

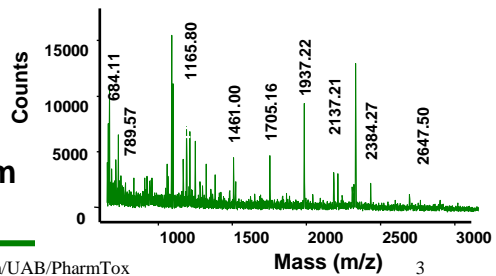

Feb 1, 2008

HelenKim/UAB/PharmTox

Types of high-throughput separation & analysis technologies

- I. 2-dimensional electrophoresis (2-DE)**
 - A. “regular” IEF/SDS-PAGE
 - B. 2D-blue-native electrophoresis
 - C. Visualization methods
- II. 2-dimensional liquid-based LC/LC**
- III. Free-flow electrophoresis**
- IV. “Chip” technology: arrays of ligands for proteins**

Feb 1, 2008

HelenKim/UAB/PharmTox

4

Parameters that govern the choice of protein separation method

- Purity of protein
- Speed of purification
- Quantity of protein
- What is the question: most important
 - Discovering a new protein/proteome
 - Identifying protein-protein interactions
 - Identifying potential modifications of known proteins

Feb 1, 2008

HelenKim/UAB/PharmTox

5

The Elements of any 2-D separation in proteomics


- The experiment!! that generates the “signal”
- Subfractionation to enrich for suspected proteins
- (Trypsin-digestion to generate peptides of the parent proteins)
- Sample work-up
- 2-D separation
- Image or other analysis to identify gel “spot” differences between untreated & treated
- MS of (trypsin-digests) spots/proteins/peptides of interest, to identify and characterize the protein

Feb 1, 2008

HelenKim/UAB/PharmTox

6

Sample preparation for 2DE:

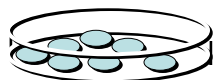

Harvest, rinse, and pellet the cells;



or



Dissect out tissue, organ, or fluids;



• Homogenize/lyse in buffer that *unfolds the proteins w/o adding or disturbing the charges*:








- High urea usually 5-8 M---unfolds the protein
- Sometimes 2 M thiourea---unfolds the protein
- 1-4% detergent---solubilizes hydrophobic components
- Beta-mercaptoethanol or other reductant, such as TBP, DTT
- Inhibitors: of proteases, kinases, & phosphatases
- Clarify by centrifugation to get rid of insoluble matter;
- Protein assay to know how much and how concentrated

Feb 1, 2008

HelenKim/UAB/PharmTox

7

Structures important in 2D electrophoresis to know and understand

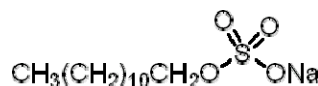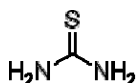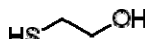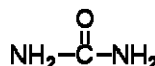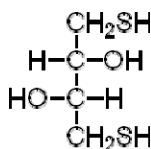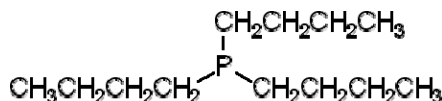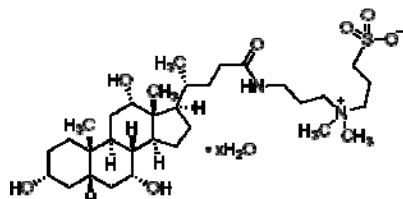


Feb 1, 2008

HelenKim/UAB/PharmTox


8

What 2-D electrophoresis involves:

- 1st dimension: Isoelectric focusing
(separation according to charge)

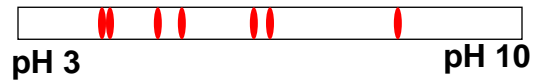

- 2nd dimension: (SDS)-PAGE
(separation according to size)

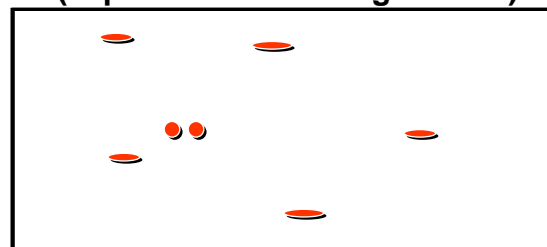

High m.w.




Low m.w.

Feb 1, 2008

HelenKim/UAB/PharmTox

9

A real 2-D gel

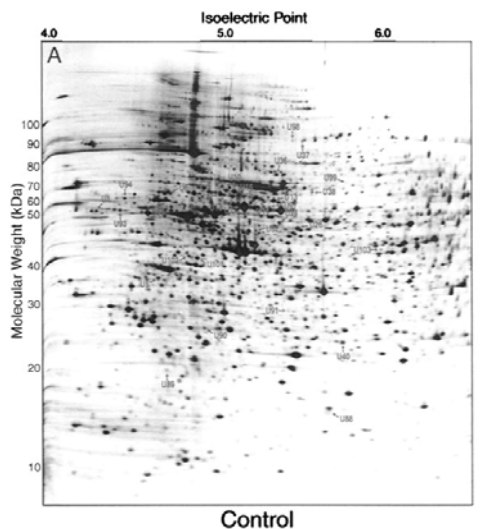

Find this and other 2-D gels at <http://www.expasy.org>

Lewis et al., [2000] Molec. Cell, 6)

(from Natalie Ahn's lab)


The pattern itself is information; a change in intensity or position of a spot has biological meaning.....similar to astronomy.

Feb 1, 2008

HelenKim/UAB/PharmTox

10

Critical part of 2-D gel proteomics: Image analysis

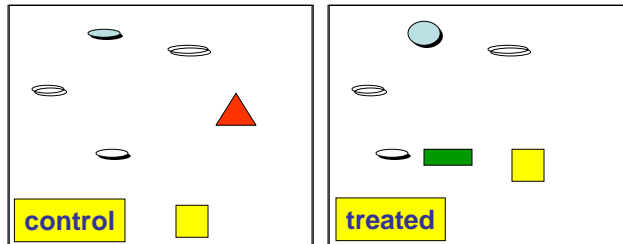

Either manually
or with software:
“compare”
the images.

Types of information:

-  ----- Upregulation of gene
-  ----- Posttranslational modification
-  ----- Downregulation of gene
-  ----- Aberrant processing

Feb 1, 2008

HelenKim/UAB/PharmTox

11

Elements of image analysis of “regular” 2D gels:

1. Compare the 2D displays of spots
2. Determine total spot number for each display
3. Quantify spot intensities, identify differences
4. Identify spots that may have “ moved” horizontally; these are candidates for those that are altered in charge, reflecting posttranslational modifications.

Ultimate and simple goal of image analysis


_____ to answer the question,
“What is changing, and by how much?”

Feb 1, 2008


HelenKim/UAB/PharmTox

12


A pair of 2D gels representing rat brain protein changes induced by ingestion of grape seed

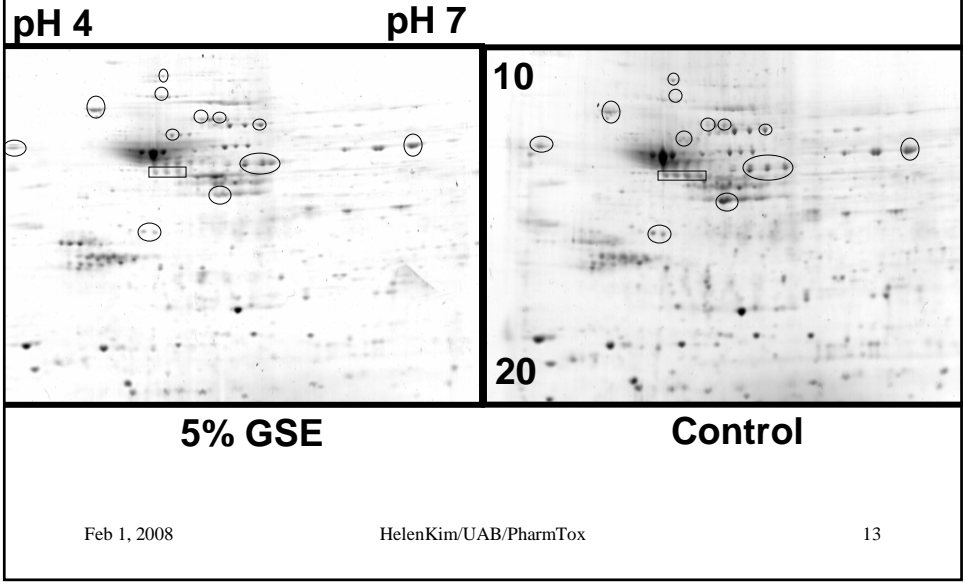

Categories of data generated by 2D gel image analysis

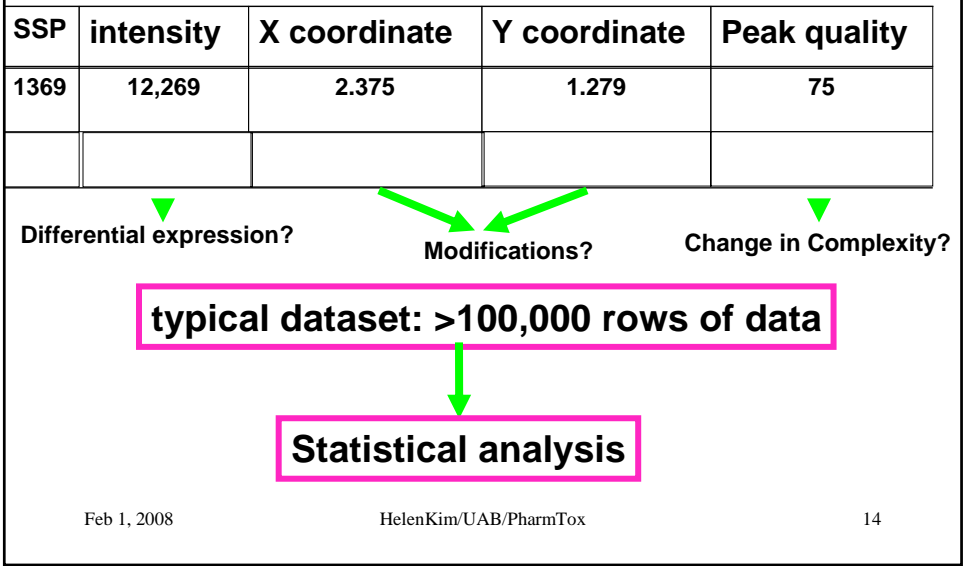

Difference gel electrophoresis (DIGE); a protein migrates with itself---- can be really valuable.

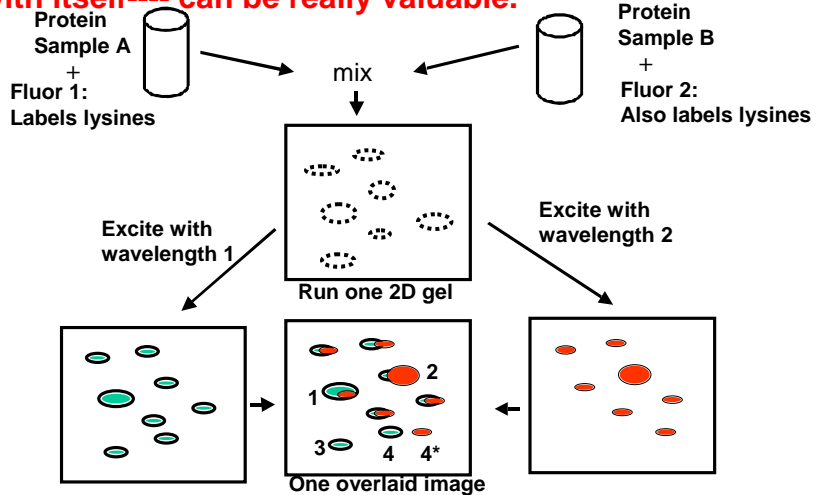

Single biggest bottleneck, even with DIGE: even commercially made gels (1st dimension, 2nd dimension) are not perfect. So, the simpler you can keep the system, the better. If you don't HAVE to use a gradient gel, don't.


Feb 1, 2008

HelenKim/UAB/PharmTox

15

Example of DIGE:

**Differential protein labeling with Cy3 and Cy5
Superimposed images from the same gel
of normal and cancer cell lines from the breast**

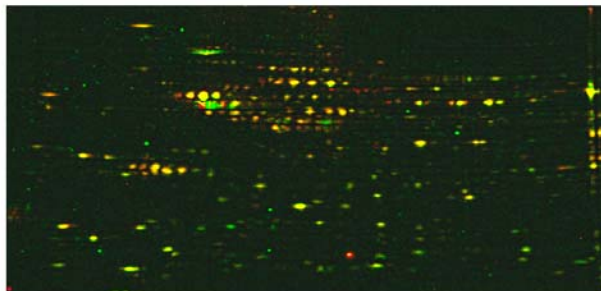

Visually: 100% green spots are specific to normal cells, 100% red spots are specific to cancer cells, degrees of yellow-orange indicate differential expression.

Feb 1, 2008

HelenKim/UAB/PharmTox

16

Typical Database obtained from proteomics analysis of breast cancer cell lines

Sample	Spot#	Protein Name	MOWSE	Accession#	Obs: kDa	Pred: kDa	Obs: pI	Pred: pI
MCF-7	BR-8	Cytoskeletal keratin-7	132	gi12803727	20.7	51.3	8.6	5.3
	BR-9	Alpha-1 anti-trypsin	90	gi1942953	19.7	44.3	9.6	
	BR-1	Human	94	gi87303	56.2	53.5	5.9	5.6
MCF-10AT		Cytoskeletal keratin-8						
	BR-6	Human	74	gi18573275	23.4	24.1	6.3	
	BR-5	hypo:XP109048 IgG heavy chain	60	gi5106591	18.1	6.4	5.7	
	BR-7	variable region	105	gi1942953	37.0	44.3	6.9	
	BR-4	Intact recomb: alpha-1 antitrypsin mutant F-L IgG heavy chain var:reg	65	gi5106591	10.0	6.4	4.6	
MCF-10A	BR-3	Human cytokeatin-8	65	gi87303	56.0	53.5	5.8	5.6
	BR-2	DNA replication silencing factor MCM-4	58	gi1705520	22.8	96.5	6	6.6

(Deshane, Johanning, and Kim, unpublished data)

Feb 1, 2008

HelenKim/UAB/PharmTox

17

“Mine” your proteomic data


- Note every difference--eventually it all means something;
- But make sure the difference is “real.”
 - What is the variation in that parameter (mw, pI) for that samespot in that treatment group;
 - Quality control issues come into play here;
 - Did you design the experiment with a statistician?
- Make sure your “basal” mw and pI are consistent with predicted and/or what others have observed;
- *Then* you can conclude that a difference in pI, for example, indicates a change in modification
- If some/all of a spot is found at a pI different from predicted, it may be constitutively modified in the “unstimulated”/“normal” group

Feb 1, 2008

HelenKim/UAB/PharmTox

18

**Metabolic labelling can enhance 2D gel analysis:
i.e. ³⁵S-methionine-labelling**

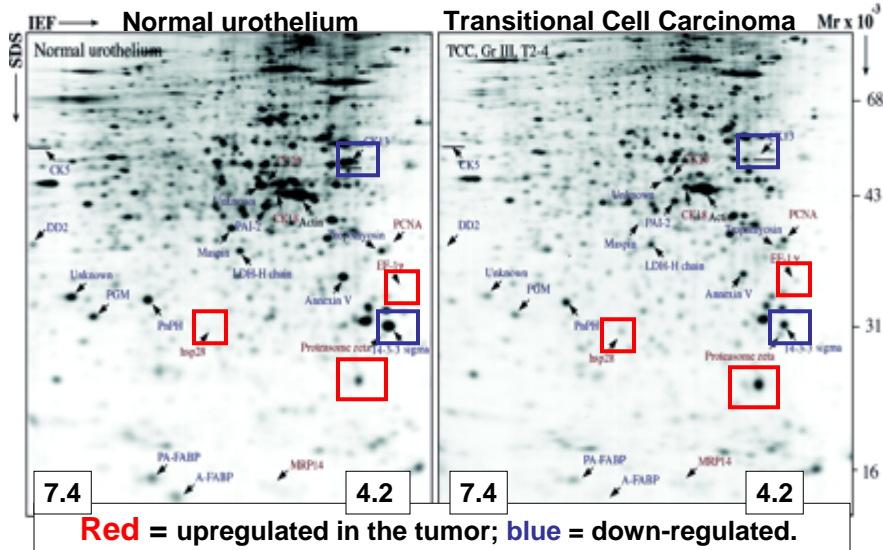

Feb 1, 2008

HelenKim/UAB/PharmTox

19

(Celis et al., Mol.Cell Proteomics, 2002, 1:4)

Visualizing 2D gels: Coomassie Brilliant Blue

QuickTime™ and a TIFF (LZW) decompressor are needed to see this picture.

In acidic conditions, the anion of CBBR combines with the protonated amino groups on proteins via electrostatic interactions.


- Inexpensive
- Image readily acquired by scanning at visible wavelengths
- No covalent mass change

Feb 1, 2008

HelenKim/UAB/PharmTox

20

Silver stain

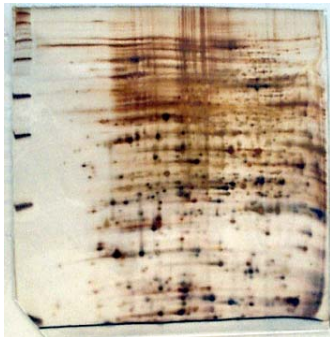

Silver ions (from silver nitrate) are chemically reduced to metallic silver on lysine residues.

This is the most sensitive protein stain, but also the least useful for quantitation because of its low dynamic range.

2D gel image courtesy of the U. Va. Mass Spectrometry Shared Facility, 2006

Feb 1, 2008

HelenKim/UAB/PharmTox



21

Sypro Ruby, and Ru(II)-BPDS Dye

Courtesy of Dr. Jim Mobley

Sypro Ruby (*fluorescent*)

- Equal to silver stain in sensitivity.
- But ~200-500-fold greater dynamic range.
- Compatible for MS analysis.

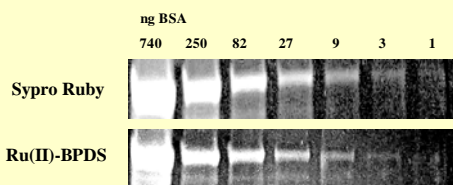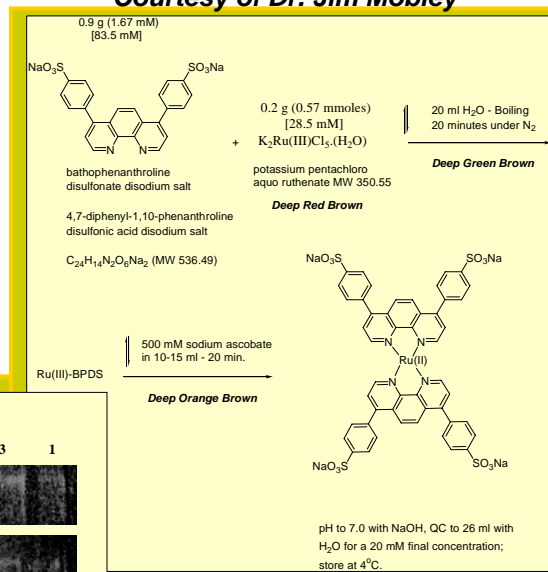

Advantages and limitations of the types of stains


	Sensitivity	Dynamic range	MS- compatible
CBBR	8 ng	10-30 x	yes
Silver	1 ng	< 10 x	Not without special precautions
Fluorescent	2 ng	3 orders of magnitude	yes

Feb 1, 2008

HelenKim/UAB/PharmTox

23

Proteins, proteins everywhere, but where's my receptor?

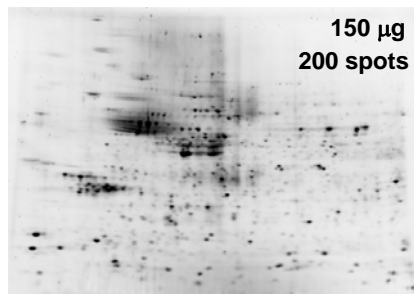

2D gel of rat brain, stained with Sypro Ruby

The genome predicts:
20,000-50,000 polypeptides.

So, 200 spots is <1% of the total proteome.

Conclusion: Even a fluorescently stained 2D gel of an unfractionated sample, only allows detection of the "low hanging fruit."

Feb 1, 2008

HelenKim/UAB/PharmTox

24

Central issue in proteome analysis: dynamic range

In cells: protein amounts vary over a wide *dynamic range*:

In blood, albumin is 3.5 g/100 ml (35 g/L = 0.5 mM) (10^{-3} M), whereas cytokines are pM (10^{-12} M)

- This is a difference of nine (9) orders of magnitude.
- A 2D gel that is overloaded with respect to an abundant protein, may have *barely detectable* amounts of a low abundance protein.....
 - If you can't see it, you don't know a protein is there;
 - Even if you know it's there, you can't do MS, because there isn't enough protein.
- **No one stain will detect 9 orders of magnitude differences in abundance of proteins.**

Feb 1, 2008


HelenKim/UAB/PharmTox

25


For greatest sensitivity, and the most biological information:

Stain the gel.....


----BUT do a Western blot or two,
with replicate gels:

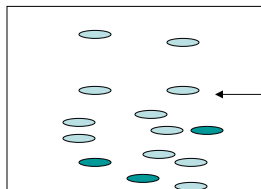

Western blot for phosphoproteins

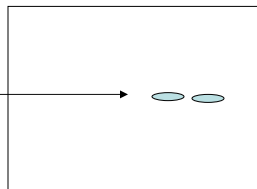

Western blot for a particular protein

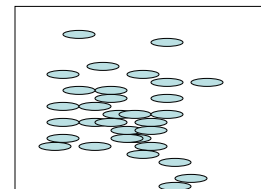

Stained gel

Feb 1, 2008

HelenKim/UAB/PharmTox

26

Deep Purple: the new fluorescence dye from GE

250

BLUE: spots on Deep Purple-stained gels;

80

RED: spots on Sypro Ruby-stained gels;

QuickTime™ and a TIFF (LZW) decompressor are needed to see this picture.

25

GREEN: spots matched to every gel

10

Deep Purple MAY detect more smaller proteins (<25kda).

Feb 1, 2008

HelenKim/UAB/PharmTox

(adapted from Duncan Veal, 2006)²⁷

Multiplex Proteomics: ProQ Emerald followed by Sypro Ruby can identify multiple glycosylated proteins at once

ProQ Emerald

Sypro Ruby

QuickTime™ and a TIFF (LZW) decompressor are needed to see this picture.



Feb 1, 2008

A- A B- B C- C

HelenKim/UAB/PharmTox



A- A B- B C- C 28

(Duncan Veal, Fluorotechnics Pty, Ltd.)

Multiplexing: valuable when sample is scarce or difficult to obtain

Normal
liver

Liver
tumor

QuickTime™ and a
TIFF (LZW) decompressor
are needed to see this picture.

Sypro Ruby:

Pro Q
Emerald:
glycosylated

Pro Q Diamond

The value of fluorescent
dyes: each can be detected
separately from others on
the gel, due to different
excitation and emission
spectra.

Feb 1, 2008

Modified from Duncan Veal,
Fluorotechnics Pty, Ltd.

SUMMARY of 2D gel stains


- **Protein stains differ according to**
 - Sensitivity/Dynamic range/MS-compatibility/Ease of capture of information
- **The fluorescent dyes offer unparalleled protein analytical capabilities due to the wide dynamic range, and their MS compatibility;**
- **“Multiplexing” allows analysis of subproteomes in the same gel, maximizing use of scarce samples;**
- **Yet, each stain has utility depending on experimental goals.**

Feb 1, 2008

HelenKim/UAB/PharmTox

30

2D-Blue-Native gels: for hydrophobic proteins

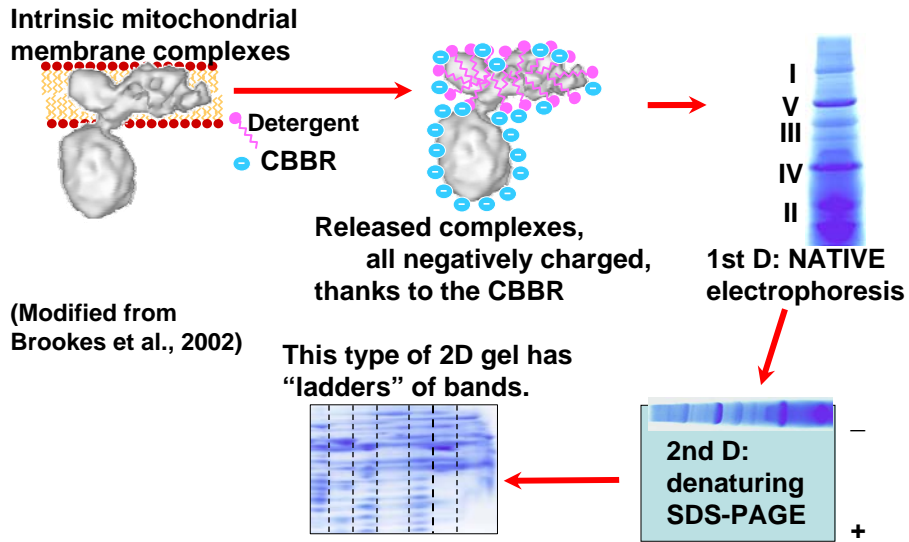


Feb 1, 2008

HelenKim/UAB/PharmTox

31

What questions does BN electrophoresis address:


(a) Which proteins are actually interacting with which?

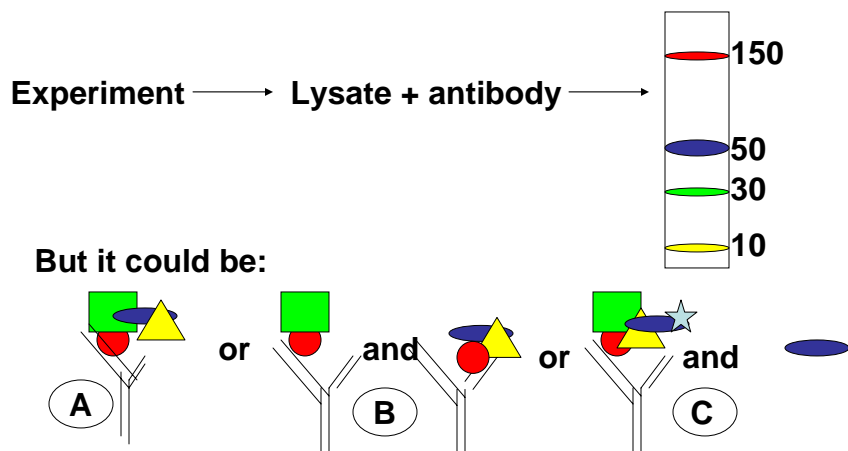



Feb 1, 2008

HelenKim/UAB/PharmTox

32

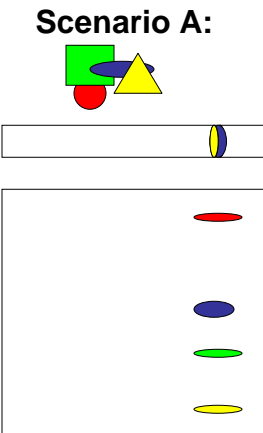Scenario A: 

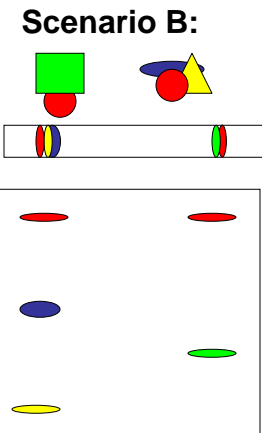Scenario B: 

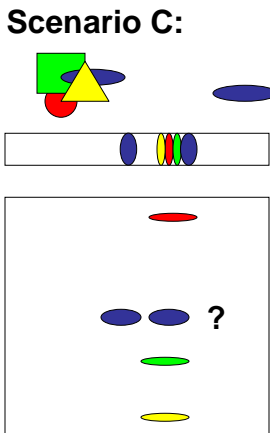
Scenario C: 

When would you see the same protein in two lanes on the 2nd?

Feb 1, 2008 HelenKim/UAB/PharmTox 33

Issues in 2D gel analysis:

I. Now you see it, now you don't:

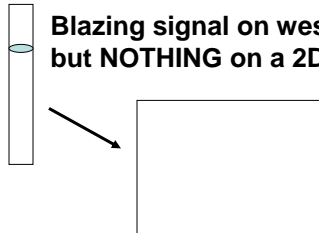 **Blazing signal on western blot of 1D SDS gel; but NOTHING on a 2D gel of the same Ip:**

Explanations?

1. Solubilization
2. Quantity
3. Modification


II. I see the band on my 8 cm long 1D gel that runs with my western blot band;
Why do I need to run a 2D gel anyway?

1. Information
2. Separate multiple proteins in that one band
- 3.

Feb 1, 2008 HelenKim/UAB/PharmTox 34

2D-LC-LC

1st dimension: chromatofocusing (like IEF, but *in solution*, so can take higher protein loads)

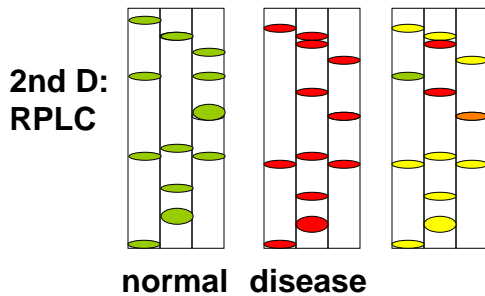


Digitized super-imposed images;
Image analysis fundamentally different from 2D gels;
Quantitation MUCH Faster;
CAVEATS?
cost;
relatively new;
reproducibility

Feb 1, 2008

HelenKim/UAB/PharmTox

35

Free flow electrophoresis

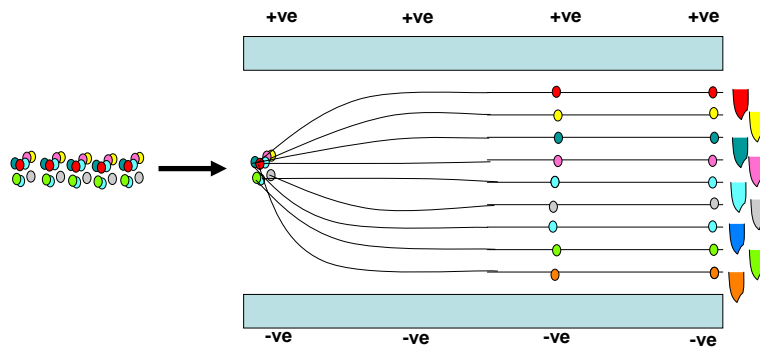

There are subpopulations of mitochondria in most cells, each with a different net charge. These can be passed through a “chamber” which has a potential difference across the sides of the chamber. Each particle will find a position where its surface charge corresponds to that of the gradient, and move along that position parallel with others moving at their own positions. This allows separation of multiple subgroups of mitochondria and physical collection into 96-well plates at the end of the free-flow chamber.

Feb 1, 2008

HelenKim/UAB/PharmTox

36

Bottom line for effective 2D gel separation and analysis


- Reduce proteome complexity by incorporating biological information or properties :
 - Intracellular location--subcellular fractions
 - Protein-protein interactions--immunoprecipitations, BN gels
 - Different states of oligomerization in vivo: microtubules

Feb 1, 2008

HelenKim/UAB/PharmTox

37

Various ways to reduce proteome complexity: Subcellular fractionation by differential centrifugation

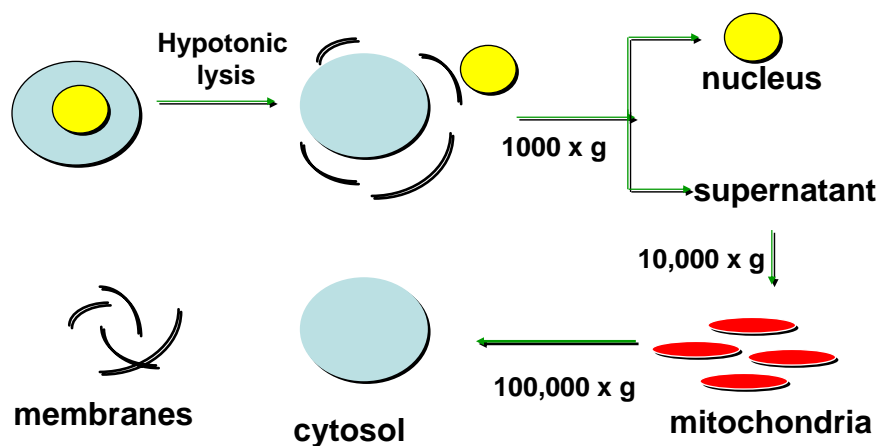

Feb 1, 2008

HelenKim/UAB/PharmTox

38

The good news: subcellular proteomes are readily “catalogued.”

Compartment	# polypeptides in SWISSPROT as of 2000
Mitochondria (1000/cell)	269-----5% of total
Lysosome (400/cell)	50-----1% of total
Peroxisome	35-----0.6%
ER and Golgi apparatus	157-----3%
Nuclei (5% cell volume)	964-----17%
Others (cytosol, membrane)	4228----75%
	total:5703


(Jung et al. [2000] Electrophoresis)
 Note date of article: this is already old data

Feb 1, 2008


HelenKim/UAB/PharmTox

39

Biological specificity of antibodies is invaluable in reducing the complexity of the proteome to be analyzed

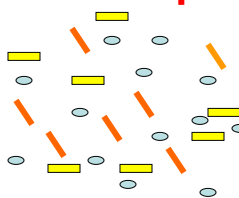

A cell lysate: 6,739 polypeptides

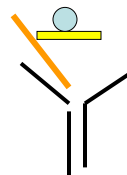

An immune complex of 1 - 3 polypeptides
 (why might there be more than one polypeptide?)

Which sample would you rather resolve on a 2D gel?

Make use of databases and the internet:

I. Check existing databases and web-links:

www.expasy.org

many are annotated

helpful links:

proteomics tools

II. Keep up with the literature/ competition:

Electrophoresis

Proteomics

Molecular & Cellular Proteomics

J. Proteome Research

III. Use genomics information when available:

The polypeptide sequence (from the cDNA) can predict electrophoretic parameters-- m.w. & pI;

helpful in setting up 2D gel conditions

Feb 1, 2008

HelenKim/UAB/PharmTox

41

Take home message

- Use of proteomics technologies enables *global* analysis of protein changes ;
- 2D electrophoresis can indicate both differential expression or posttranslational modifications;
- Choice of separation governed by
 - Abundance of sample
 - Question being asked
 - Technology available to you
 - Cost is a factor;

Feb 1, 2008

HelenKim/UAB/PharmTox

42

Future directions in intact protein analysis approaches

- I. DIGE and Cy-dye labelling will enhance 2D gel analysis of complex proteomes;**
- II. Subcellular fractionation will regain importance in proteome analysis;**
- III. While automated 2D LC/LC-MS/MS may appear more highthroughput for “discovery,” every resolved spot on a 2D gel is a discovery;**
- IV. 2D gel positional information, *without protein identities*, is information itself.**
- V. Where 2D gels may play larger roles is in validation of results generated by other approaches;**

Feb 1, 2008

HelenKim/UAB/PharmTox

43

Final thoughts

- What proteomics technology gives back is like any other analytical approach: it's as good as what you put in;**
- Be mindful of distinguishing between low abundance proteins vs low level contaminants;**
- Keep in mind “conventional” approaches like Western blotting to validate proteomic results;**
- Purify, purify, purify before running any proteomic experiment.**
- Formulating a hypothesis forces you to incorporate biological information that can enhance proteomic analysis.**

Feb 1, 2008

HelenKim/UAB/PharmTox

44